

Sample: 240610-5Mo1.5FeACcal-IWI-mk-021
Operator: Bryce
Submitter:
File: D:\Flex\data\Bryce\000-051 240610 5Mo1.5FeAC_MK021 H2 TPR.SMP

Started: 6/10/2024 11:46:32 AM
Completed: 6/10/2024 4:21:08 PM

Sample mass: 0.1080 g
Report time: 6/12/2024 4:08:33 PM

Comments: Cold trap in LN2/IPA slurry
TPR with 10% H2
5Mo1.5FeAC IWI MK021

Summary Report

Experiment 1: AgO TPR with H2Ar

Analysis type: Temperature Programmed

Calibration: None

Measured flow rate: 2,230.8 $\mu\text{mol}/\text{min}$

Signal offset: 0.00000

Signal inverted: No

Peak Number	Temperature at Maximum (°C)	Area	Peak Height
1	201.6	1.398e-03	2.242e-04
2	488.7	0.95791	-6.900e-02


Sample: 240610-5Mo1.5FeACcal-IWI-mk-021
Operator: Bryce
Submitter:
File: D:\Flex\data\Bryce\000-051 240610 5Mo1.5FeAC_MK021 H2 TPR.SMP

Started: 6/10/2024 11:46:32 AM
Completed: 6/10/2024 4:21:08 PM

Sample mass: 0.1080 g
Report time: 6/12/2024 4:08:33 PM

Comments: Cold trap in LN2/IPA slurry
TPR with 10% H2
5Mo1.5FeAC IWI MK021

Signal (a.u.) vs. Time

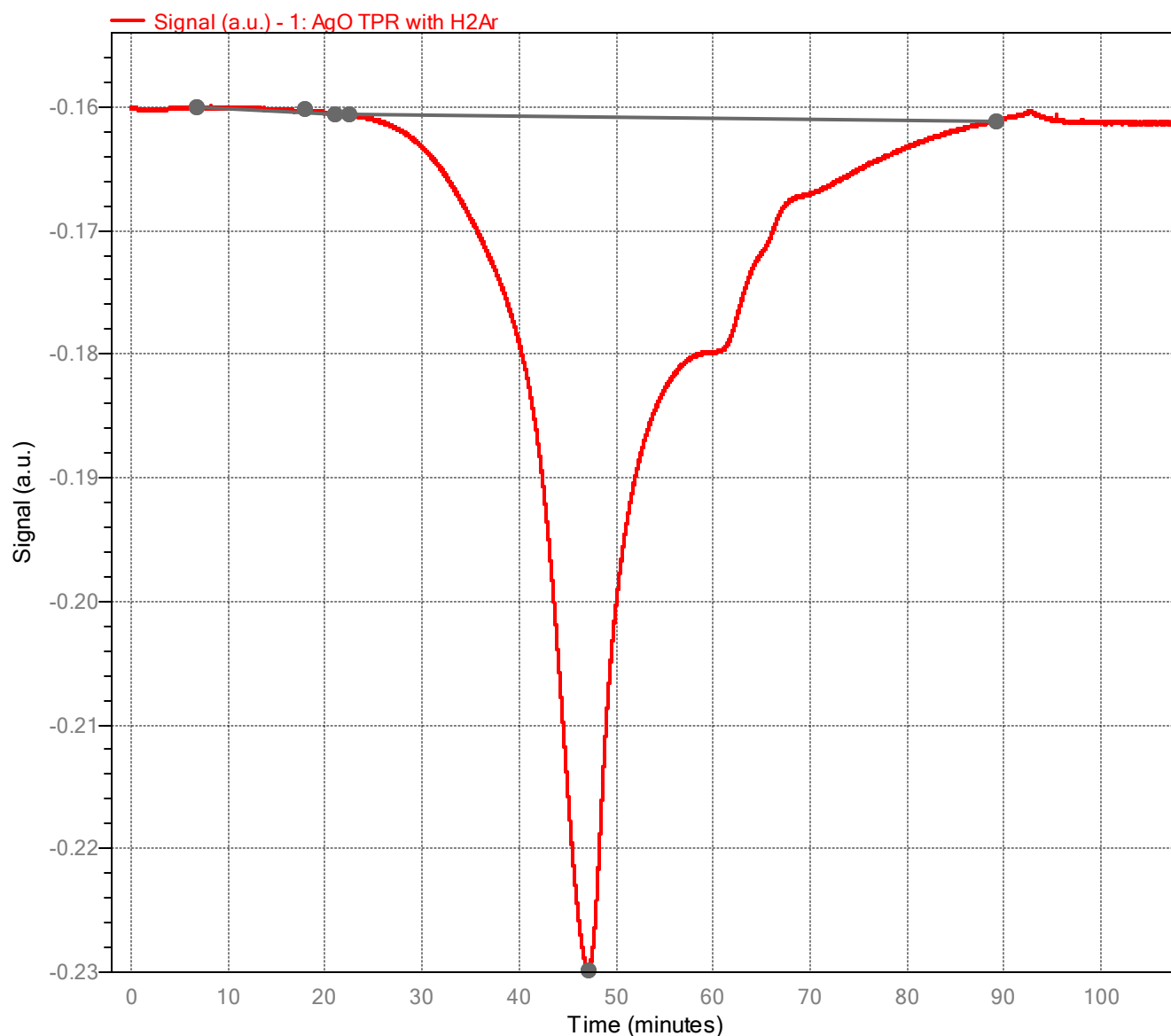


Sample: 240610-5Mo1.5FeACcal-IWI-mk-021
Operator: Bryce
Submitter:
File: D:\Flex\data\Bryce\000-051 240610 5Mo1.5FeAC_MK021 H2 TPR.SMP

Started: 6/10/2024 11:46:32 AM
Completed: 6/10/2024 4:21:08 PM

Sample mass: 0.1080 g
Report time: 6/12/2024 4:08:33 PM

Comments: Cold trap in LN2/IPA slurry
TPR with 10% H2
5Mo1.5FeAC IWI MK021

Temperature vs. Time

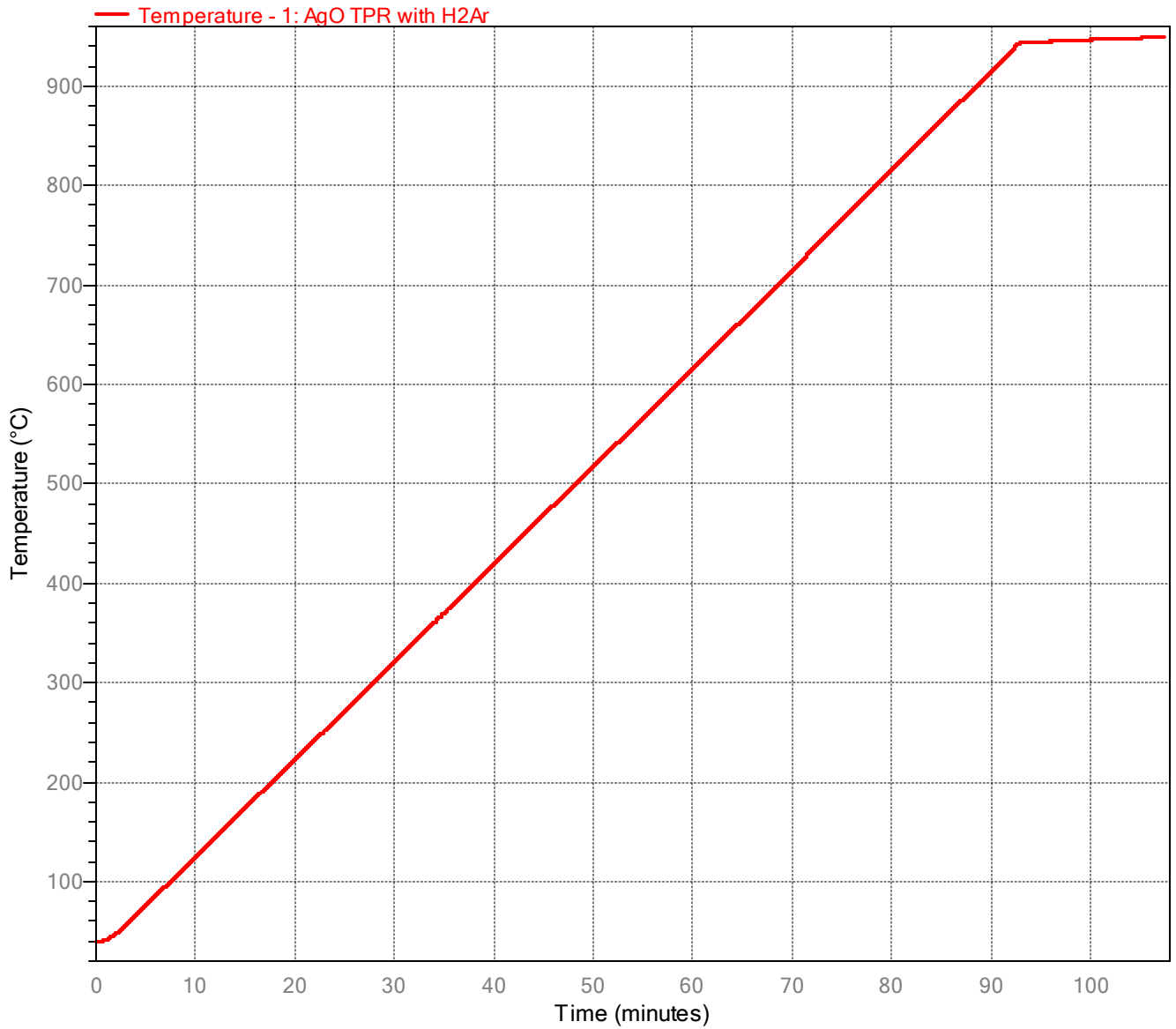


Sample: 240610-5Mo1.5FeACcal-IWI-mk-021
Operator: Bryce
Submitter:
File: D:\Flex\data\Bryce\000-051 240610 5Mo1.5FeAC_MK021 H2 TPR.SMP

Started: 6/10/2024 11:46:32 AM
Completed: 6/10/2024 4:21:08 PM

Sample mass: 0.1080 g
Report time: 6/12/2024 4:08:33 PM

Comments: Cold trap in LN2/IPA slurry
TPR with 10% H2
5Mo1.5FeAC IWI MK021

Signal (a.u.) vs. Temperature

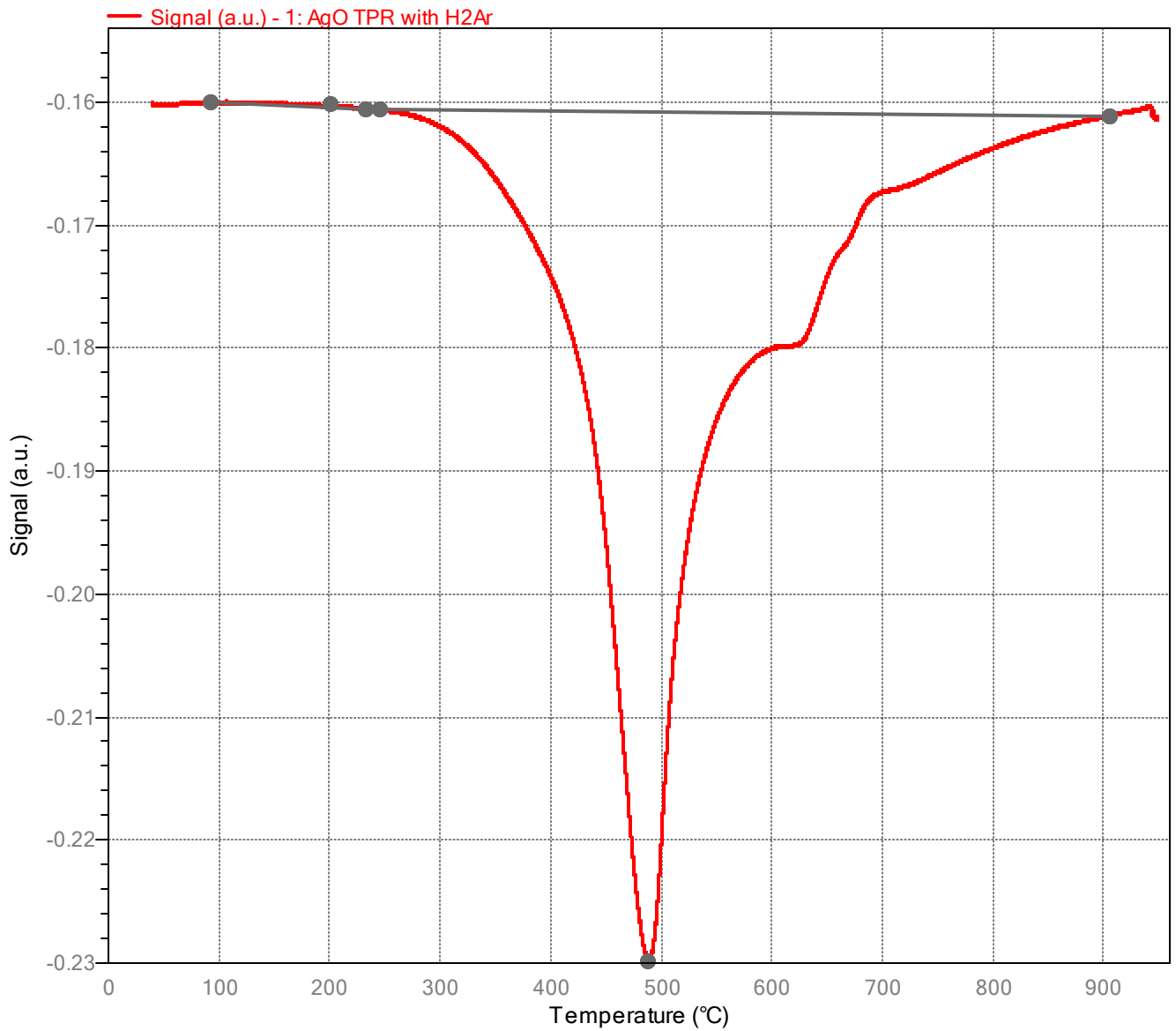

Sample: 240610-5Mo1.5FeACcal-IWI-mk-021
 Operator: Bryce
 Submitter:
 File: D:\Flex\data\Bryce\000-051 240610 5Mo1.5FeAC_MK021 H2 TPR.SMP

Started: 6/10/2024 11:46:32 AM Sample mass: 0.1080 g
 Completed: 6/10/2024 4:21:08 PM Report time: 6/12/2024 4:08:33 PM

Comments: Cold trap in LN2/IPA slurry
 TPR with 10% H2
 5Mo1.5FeAC IWI MK021

Sample Information

Method: \$ AgO TPR with H2Ar
 Sample: 240610-5Mo1.5FeACcal-IWI-mk-021
 Operator: Bryce
 Submitter:
 Mass type: Calculated
 Empty tube: 32.3797 g
 Sample + tube: 32.4877 g
 Sample mass: 0.1080 g
 Density: 1.000 g/cm³
 Type of data: Automatically collected
 Instrument type: 3500 Dynamic Chemi
 Original instrument type: 3500 Dynamic Chemi
 Comments: Cold trap in LN2/IPA slurry
 TPR with 10% H2
 5Mo1.5FeAC IWI MK021

Element	Atomic Weight	Atomic Cross Sect. Area (nm ²)	Active Metals Table				MxOy X	MxOy Y
			Density (g/cm ³)	Percent of Sample Mass (%)	Percent Reduced (%)			
molybdenum	95.940	0.0730	10.220	5.000	100.00	1	0	
	Adsorptive Stoichiometry							
	Hydrogen	2.000						
	Oxygen	2.000						
	Carbon Monoxide	1.000						
Iron	55.847	0.0613	7.890	1.500	100.00	1	0	

Analysis Conditions

Analysis conditions: 10% H2
 View conditions for:

Baseline Options
 Stable Baseline
 Slope threshold: 0.010 %/min
 Duration: 5.00 min
 Change from Baseline
 Acceleration threshold: 0.200 %/min²
 Duration: 0.10 min
 Return to Baseline
 Acceleration threshold: 0.050 %/min²

Sample: 240610-5Mo1.5FeACcal-IWI-mk-021
 Operator: Bryce
 Submitter:
 File: D:\Flex\data\Bryce\000-051 240610 5Mo1.5FeAC_MK021 H2 TPR.SMP

Started: 6/10/2024 11:46:32 AM
 Completed: 6/10/2024 4:21:08 PM

Sample mass: 0.1080 g
 Report time: 6/12/2024 4:08:33 PM

Comments: Cold trap in LN2/IPA slurry
 TPR with 10% H2
 5Mo1.5FeAC IWI MK021

Baseline Options

Duration: 1.00 min

Tasks

Description	Tasks	Details
1.01 Experiment	AgO TPR with H2Ar	
	Type: Temperature Programmed	
1.02 Gas Flow	Carrier gas: Nitrogen	
	Loop or injection gas:	
1.03 Temperature Ramp	Type: Sample	
	Temperature: 150.0 °C	
	Ramp rate: 10.0 °C/min	
	Hold time: 45.00 min	
1.04 Temperature Ramp	Type: Sample	
	Temperature: 40.0 °C	
	Ramp rate: 10.0 °C/min	
	Hold time: 30.00 min	
1.05 Gas Flow	Carrier gas: Hydrogen-Argon	
	Loop or injection gas:	
1.06 Set Detector	Enable detector: Yes	
	Filament temperature: 175.0 °C	
1.07 Wait	Wait until baseline is stable.	
1.08 Start Recording	One measurement every 0.10 s	
1.09 Temperature Ramp	Type: Sample	
	Temperature: 950.0 °C	
	Ramp rate: 10.0 °C/min	
	Hold time: 0.00 min	
1.10 Stop Recording Termination	Stop Recording	
	Return to ambient: Yes	
	Carrier gas:	
	Loop gas:	

Adsorptive Properties

Adsorptive: Hydrogen-Argon (H2Ar)
 Maximum manifold pressure: 925.00 mmHg
 Therm. tran. hard-sphere diameter: 0.31480 nm
 Molecular cross-sectional area: 0.143 nm²
 Adsorbate molecular weight: 39.95
 Mass flow constant: 1.000
 Relative thermal conductivity: 1.40
 Gas blend: Yes
 Inert gas: Argon
 Active gas: Hydrogen
 % Active Concentration 10.000 %
 Fluid properties: H:\FPI\fpiv6\argon.fpi
 Dosing method: Normal

Sample: 240610-5Mo1.5FeACcal-IWI-mk-021
Operator: Bryce
Submitter:
File: D:\Flex\data\Bryce\000-051 240610 5Mo1.5FeAC_MK021 H2 TPR.SMP

Started: 6/10/2024 11:46:32 AM
Completed: 6/10/2024 4:21:08 PM

Sample mass: 0.1080 g
Report time: 6/12/2024 4:08:33 PM

Comments: Cold trap in LN2/IPA slurry
TPR with 10% H2
5Mo1.5FeAC IWI MK021

Adsorptive: Nitrogen (N2)
Maximum manifold pressure: 925.00 mmHg
Therm. tran. hard-sphere diameter: 0.36810 nm
Molecular cross-sectional area: 0.162 nm²
Adsorbate molecular weight: 28.01
Mass flow constant: 1.000
Relative thermal conductivity: 1.00
Gas blend: No
Fluid properties: H:\FPI\fpiv6\nitrogen.fpi
Dosing method: Normal

Report Options

Peak Detection/Integration Options
Baseline mode: Best fit baseline
Minimum peak height: 0.25 % F.S.
Peak smoothing: 0 points
Sensitivity: 1.0e-02 % F.S.
Maximum shoulder ratio: 33%
Maximum group separation: 33%
Minimum peak area: 1.0 % F.S.·min
Maximum baseline slope: 0.1 % F.S./sec

Sample: 240610-5Mo1.5FeACcal-IWI-mk-021
 Operator: Bryce
 Submitter:
 File: D:\Flex\data\Bryce\000-051 240610 5Mo1.5FeAC_MK021 H2 TPR.SMP

Started: 6/10/2024 11:46:32 AM
 Completed: 6/10/2024 4:21:08 PM

Sample mass: 0.1080 g
 Report time: 6/12/2024 4:08:33 PM

Comments: Cold trap in LN2/IPA slurry
 TPR with 10% H2
 5Mo1.5FeAC IWI MK021

Sample Log

Date	Time	Log Message
6/10/2024	11:46:32 AM	1.01 Experiment started. Type: Temperature Programmed
6/10/2024	11:46:32 AM	1.02 Set gas flow. Carrier gas: N2 Flow rate: 50.00 cm ³ STP/min Loop gas:
6/10/2024	11:47:22 AM	1.03 Ramping Sample temperature at 10.0 °C/min to 150.0 °C for 45 minutes.
6/10/2024	12:48:58 PM	1.04 Ramping Sample temperature at 10.0 °C/min to 40.0 °C for 30 minutes.
6/10/2024	1:33:08 PM	1.05 Set gas flow. Carrier gas: H2Ar Flow rate: 50.00 cm ³ STP/min Loop gas:
6/10/2024	1:33:08 PM	1.06 Set detector. Detector enabled: Yes Filament temperature: 175.0 °C
6/10/2024	1:33:08 PM	1.07 Wait for baseline stable. Time: 300 s
6/10/2024	1:54:59 PM	1.08 Recording started. Carrier gas: H2Ar Flow rate: 50.00 cm ³ STP/min Loop gas: Recording period: 0.10 s
6/10/2024	1:54:59 PM	1.09 Ramping Sample temperature at 10.0 °C/min to 950.0 °C for 0 minutes.
6/10/2024	3:42:24 PM	1.10 Recording ended.
6/10/2024	3:42:24 PM	2.01 Analysis terminated. Detector enabled: No Return to ambient temperature: Yes
6/10/2024	4:21:08 PM	Finished a sample analysis for D:\Flex\...000-051 240610 5Mo1.5FeAC_MK021 H2 TPR.SMP.